


User Instructions

VH1003 - Over Bed Table

VH1001 - Divan Over Bed Table

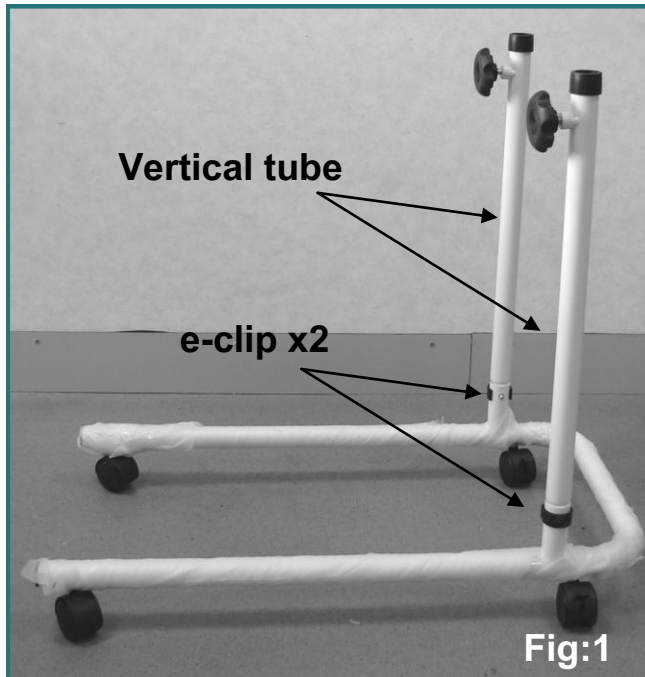

Verax Healthcare Ltd, Unit 1, Brook Mill
Industrial Estate, Wrea Green, Preston,
Lancashire, PR4 2PH

Email: sales@veraxhealthcare.com

www.veraxhealthcare.com

General Safety Warning

Please ensure these instructions are fully read and implemented. Failure to do so may result in injury to the user.

Retain in a safe place for future reference.

Intended Use

The Wheeled Overbed/ Divan Overbed Tables are intended for use when sitting in either a chair or bed, to provide an adjustable height/angle table surface for reading, eating etc. It is not to be used as a trolley or designed to be weight bearing.

Contents

Contents	Qty
Base Frame*	1
Vertical tubes with pin e-clips	2
Table top with 2 telescopic tubes	1

*VH1003 with castors (illustrated),
VH1001 without (not illustrated).

Assembly and Height Adjustment

1. Unpack all contents, check all parts are present and in good condition before attempting to assemble the table.
2. Assemble the vertical tubes into the base and secure using the pin e-clips. See Fig. 1
3. Slide the two telescopic tubes attached to the table top into the two open tube ends of the side frames (Fig. 2)
4. Adjust the height and angle of the table top as required and tighten hand wheels clockwise to secure (Fig. 2 Inset)

Safety Precautions

- ☒ Inspect at regular intervals for wear and damage.
- ☒ Do not overload table top (10Kg Max load)
- ☒ Check hand wheels and e-clips are secure at regular intervals.

Cleaning Instructions

Clean regularly, wiping down with mild detergent /disinfectant solution.

The Overbed Table is part of a wide range of daily living aids available from Verax Healthcare.